

Bibliographic Fields

(54)

考案の名称 吸尿テープ

(21)

実願 昭 59-45905

(22)出願

昭 59(1984)3 月 31 日

(72)

考案者 横山小枝子 与野市大戸 705

(71)出願人

横山小枝子 与野市大戸 705 明細書

Abstract**明細書****1. 考案の名称**

吸尿テープ

Claims**2. 実用新案登録請求の範囲**

赤ちゃん又は病人が放尿した時におしめを解いて調べなくとも、本考案の吸尿テープを性器の付近に設置しておけば排尿したかどうかは、この考案テープで尿の吸い加減でわかる。

Specification**3. 考案の詳細な説明****イ. 考案の目的**

リトマス試験紙のような吸取紙で尿を吸い上げ、テープの中の吸取紙(第1図の2が変色していくさまで排尿後か、否か知ることが出来る。

ロ. 考案の構成

図面の1の材質は、透明で薄いビニール製の平たい筒状のもの

図面の2の特質は、リトマス試験紙のような素材で吸収力の強力なもので変色して変わっていくもの。

第1図の3の部分は、吸取紙がむき出しになっているところに脱脂綿のようなもので覆うことによって尿を少し多く含ませ吸取紙の吸収効果を高める。

この本考案テープの長さは、性器からズボンの

(54)

***title *urine tape

(21)

Japan Patent Application Showa 59-45905

(22) Application

*59 (1984) March 31 *

(72)

inventor Yokoyama twig *Yono City Oto 705

(71) Applicant

Yokoyama twig *Yono City Oto 705specification

Specification**1. title of device**

Sucking/absorbing urine tape

2. Utility Model Claims

When baby or patient releasing urine doing, solving the diaper, also without inspecting, if it installed sucking/absorbing urine tape of this utility model in vicinity of characteristic vessel, whether or not anti-urine it did, with this device tape you understand with inhaling allowance of urine.

3. detailed description of device**objective of a. device**

2 of Figure 1 does discoloration, it can know whether with after the anti-urine, whether no. configuration of b. device To suck up urine with collecting paper like litmus test paper, collecting paper in the tape

As for material of 1 of drawing, with transparent flat of thin vinyl make wants thing of cylinder which

As for characteristic of 2 of drawing, with material like litmus test paper discoloration doing with strong ones of absorption power, those which change.

portion of 3 of Figure 1 urine raises absorbing effect of the implication collecting paper a little more by fact that it covers with those like absorbent cotton in place where collecting paper becomes bare.

As for length of this this utility model tape, effect of length.

ウエスト部分からちよつとのぞく程度の長さ。ハ、
考案の効果

第 1 図の表側の下の部分を少し吸取紙をむき出しにし、そこに脱脂綿のようなもので覆うことによってビニールパイプの中に通っている吸取紙の吸収を高め、体温計の水銀のように上ってくるさまがよく解る。

本考案の吸尿テープの使用について、男子は、本考案テープを体の前部に当てがい(性器とオムツの間に挟む。)

女子は、本考案テープを体の後部に当がう。(お尻に敷くように挟む。)

Drawings

4. 図面の簡単な説明

第 1 図は本考案の一実施例を示す吸尿テープの表側

第 1 図のイ、ロの点線は切断箇所。

場 2 図は同上の裏側。

穿 3 図は、第 1 図イ、ロ部分の切断図。

1 はビニールの筒状のもの。

2 はリトマス試験紙のような吸取紙。

3 は脱脂綿。

実用新案登録出願人 横山小枝子

c. device of extent just a little of excluding from waist part amount of trousers from the characteristic vessel

portion under front side of Figure 1 it designates collecting paper a little as bare, it raises absorption of collecting paper which by fact that it covers with those like absorbent cotton there it is passing in vinyl pipe, like mercury of body thermometer it rises, you understand well.

Concerning use of sucking/absorbing urine tape of this utility model, as for male, this utility model tape applying is in front portion of body, (It puts between to characteristic vessel and between diaper.)

As for female, this utility model tape in rear part of body this *. (In order to spread in rear end, it puts between.)

4. Brief Explanation of the Drawing (s)

As for Figure 1 front side of sucking/absorbing urine tape which shows one Working Example of this utility model

[i] As for dotted line of [ro] of Figure 1 cutting site.

As for place 2 figures back side. of same as above

* As for 3 figures, Figure 1 [i], cross section. of [ro] portion

As for 1 those of cylinder of vinyl.

As for 2 collecting paper. like litmus test paper


As for 3 absorbent cotton.

Utility Model applicant Yokoyama twig child


第3圖



第2圖

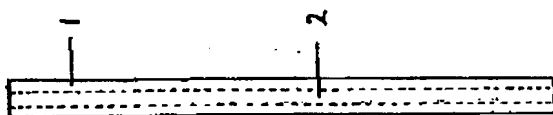

第1圖

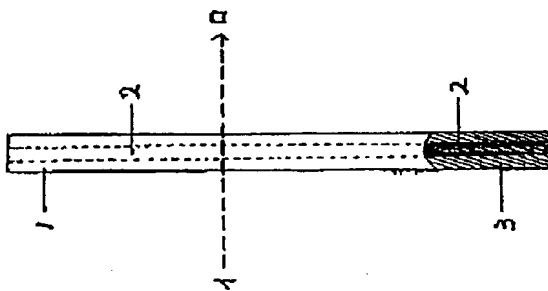

使用新案登錄人 橫山小破子

61
其附 60-161108